


Accuracy Check Label Settings

When using the **Accuracy Check** module, a QR code will print on your item and bag labels, as well as the new picking list. You may need to adjust the standard position of the QR codes if you use pre-printed labels.

You can adjust the position of the QR codes from the **Pharmacy Details - Accuracy Check**.

Set your QR Code Positions

1. From the Pharmacy Manager **Toolbar**, select **Tools - System Settings - Pharmacy Details - Accuracy Check** to display the **Accuracy Check** settings screen.

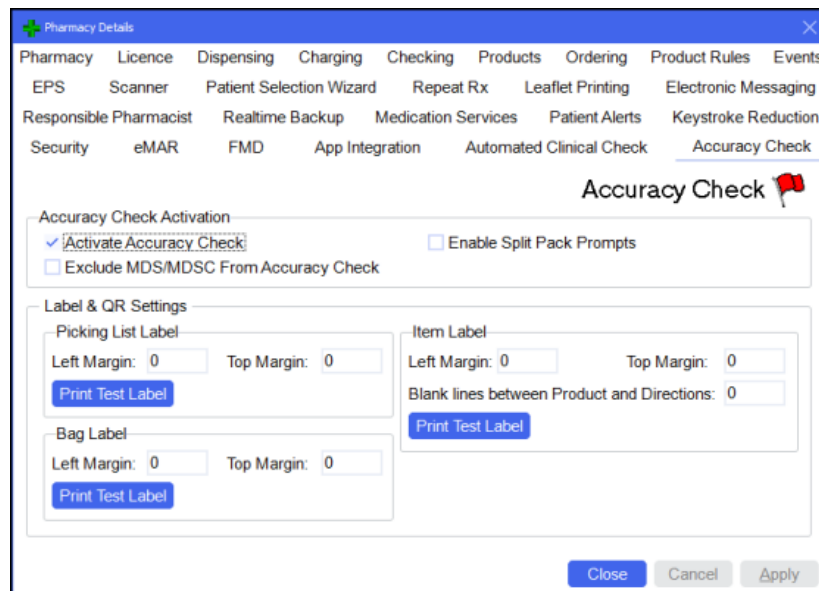

2. **Label & QR Settings** displays the current QR code positions for the picking list, item and bag labels. From here you can adjust the top and left margins to adjust the position of the code by pixels.

Please note a value of 0 for both the top and left margin defaults the QR code to the following position.

Label	Default Position of QR code when set to 0,0
Picking List	Top right
Bag	Top right
Item	Top left


Adjusting the left margin value from 0 to 10 moves the QR code to the right by 10 pixels, or adjusting the top margin from 0 to 10 moves the QR code down by 10 pixels. Alternatively, negative numbers move the QR code in the opposite direction, to the left or up.


3. Select **Print Test Label**  to print a test label for each label to check your settings.


Check your Item Labels

Due to the addition of the QR codes and using the Windows Printer Driver to print your labels, you may need to adjust your label printer configuration.


The label top margin is the gap between the top of the label and the start of the text. This needs to be changed if your text is disappearing off the top of the label or the label is pre-printed at the top and your text needs moving down.

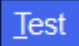1. From the **Pharmacy Manager Toolbar**, select **Tools - System Settings - Printer Configuration - Pharmacy Manager**.
2. Select the **Workstations** tab, highlight the **Label Printer** and then select **Test**  to print a test item label.

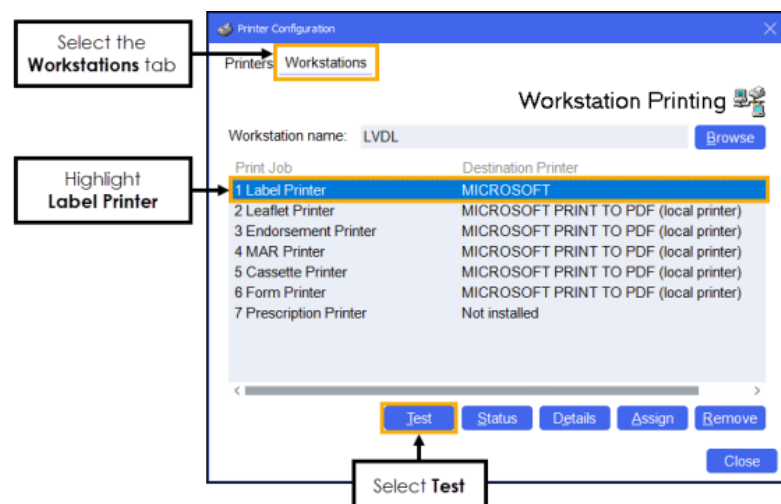

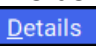3. Should you need to adjust the settings, with the label printer still highlighted, select **Details**  to open the **Printer Details** screen.

4. The **Printer Details** screen displays.

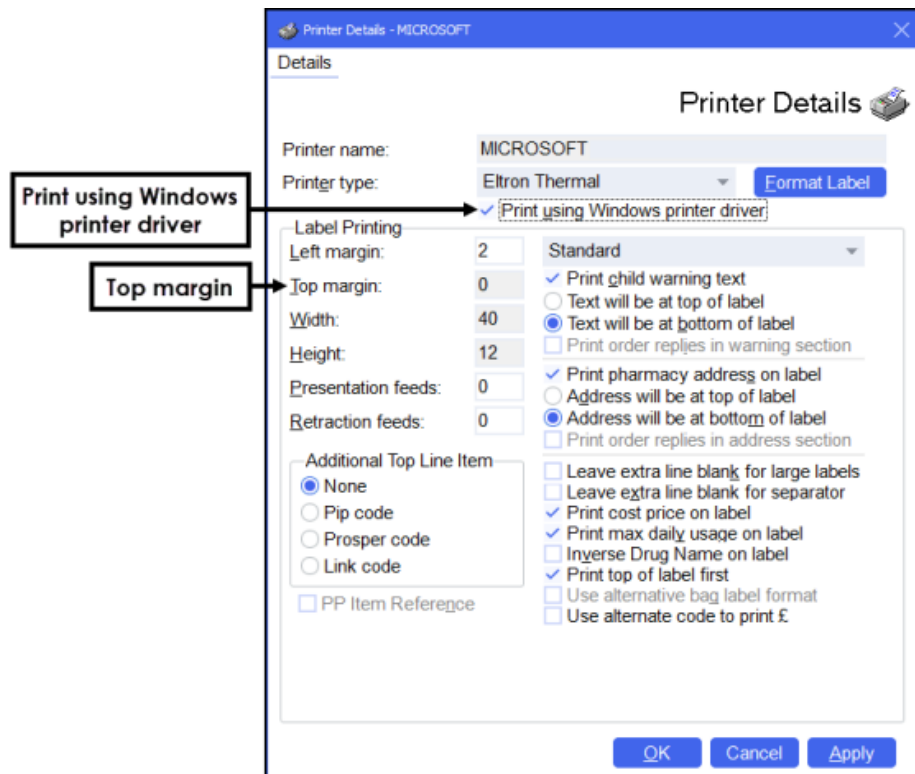

To change the top margin, **Print using Windows printer driver** must be unticked.

5. Edit **Top margin** as required. As an example, adjusting the value 0 to 1 adds a margin of 1 line of text.

6. Tick **Print using Windows printer driver** and then select **Apply** and **Close**.




i Important - No additional label formatting is possible with **Accuracy Check** enabled, therefore the label defaults to the **Pharmacy Manager** standard.

7. From the **Printer Configuration** screen, select **Test** to print another test label to check the settings and repeat as required.



Example Labels

The label examples below have used blank labels and the following settings:

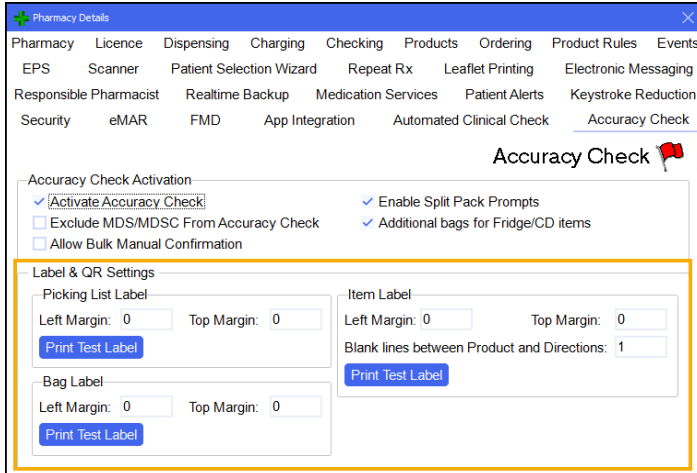

Pharmacy Details - Accuracy Check

Picking List Label

- Left Margin: 0
- Top Margin: 0

Bag Label

- Left Margin: 0
- Top Margin: 0

Item Label

- Left Margin: 0
- Top Margin: 0
- Blank lines between Product and Directions: 1

Pharmacy Details - Accuracy Check

Picking List Label

- Left Margin: 0
- Top Margin: 0

Bag Label


- Left Margin: 0
- Top Margin: 0

Item Label

- Left Margin: 0
- Top Margin: 0
- Blank lines between Product and Directions: 1

Printer Details

- Print using Windows printer driver: Ticked
- Top margin: 1

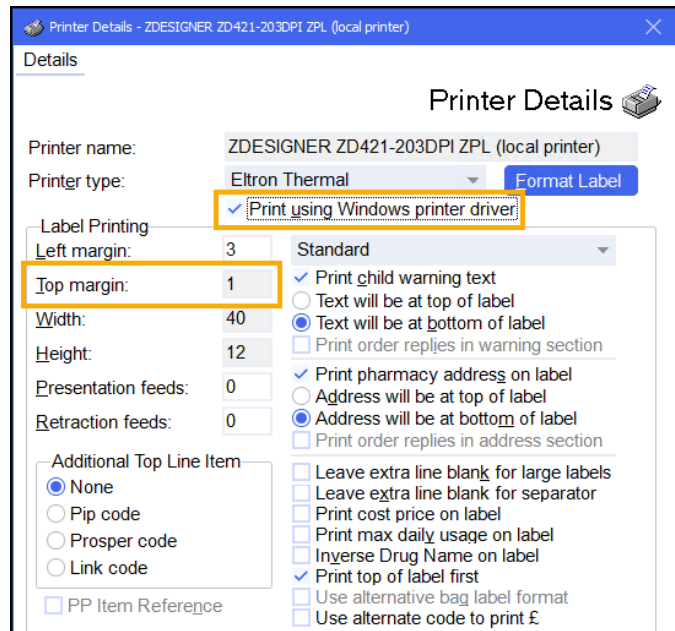

Printer Details - ZDESIGNER ZD421-203DPI ZPL (local printer)

Label Printing

- Left margin: 3
- Top margin: 1
- Width: 40
- Height: 12
- Presentation feeds: 0
- Retraction feeds: 0


Print using Windows printer driver

Additional Top Line Item

- None
- Pip code
- Prosper code
- Link code

PP Item Reference


Picking List Label

Mr Gary Green 

6215667544 (1 of 1)


ITEMS REQUIRED
24 Ibuprofen 400mg Tablets
28 Omeprazole 20mg g/r capsules

Bag Label

MR GARY GREEN
11 Short Lane
Leyland
PR26 1AP
03/02/25 2 items NHS 

KEEP OUT OF REACH AND SIGHT OF CHILDREN
Cegedim Rx Parkway 2 Buckshaw PR7 7NR

Item Label

 **24 IBUPROFEN 400MG TABLETS**

ONE to be taken **THREE** times a day when
required
Take with or just after
food, or a meal.
Mr Gary Green(3) 03/02/2025

KEEP OUT OF REACH AND SIGHT OF CHILDREN
Cegedim Rx Parkway 2 Buckshaw PR7 7NR