



# SIS 10500 Vision Release Guide

Version 2

15/04/2020

Vision  
The Bread Factory  
1a Broughton Street  
London SW8 3QJ



Registered No: 1788577 England

[www.visionhealth.co.uk](http://www.visionhealth.co.uk)  
T +44(0) 20 7501 7000  
F +44(0) 20 7501 7100

**© 2020 In Practice Systems Limited**

All Rights Reserved.

No part of this document may be photocopied, reproduced, stored in a retrieval system or transmitted in any form or by any means, whether electronic, mechanical, or otherwise, without the prior written permission of In Practice Systems Limited.

No warranty of accuracy is given concerning the contents of the information contained in this publication. To the extent permitted by law, no liability (including liability to any person by reason of negligence) will be accepted by In Practice Systems Limited, its subsidiaries or employees for any direct or indirect loss or damage caused by omissions from or inaccuracies in this document.

In Practice Systems Limited reserves, the right to change without notice the contents of this publication and the software to which it relates.

Product and company names herein may be the trademarks of their respective owners.

**In Practice Systems Limited**  
**The Bread Factory**  
**1a Broughton Street**  
**London**  
**SW8 3QJ**

Website: [www.visionhealth.co.uk](http://www.visionhealth.co.uk)



# Contents

---

<b>SIS 10500 VISION RELEASE GUIDE</b>	<b>1</b>
<b>CONTENTS</b>	<b>IV</b>
<b>SIS 10500 VISION RELEASE GUIDE</b>	<b>1</b>
<b>Summary of Changes</b>	<b>1</b>
<b>Vision+ Templates</b>	<b>3</b>
Viewing data	3
<b>Vision+ Template Designer</b>	<b>5</b>
Adding Active Repeats to a Template	5
Adding an Acute Row	7
Adding a Repeat Usage Row	8
Adding a Problem Row	10

# SIS 10500 Vision Release Guide


---

## Summary of Changes

17/03/2020

SIS 10500 delivers the following Vision+ changes:

### Vision+ SNOMED CT







-  **SNOMED CT** – The latest **SNOMED CT** dictionary is implemented in SIS 10500 and includes the Coronavirus concepts. For more information on Vision's Coronavirus guidance see - <https://info.visionhealth.co.uk/coronavirus-advice>



**Note** – All countries (excluding Northern Ireland) now have access to the Coronavirus clinical terms in Vision 3 (using local terms). However, in the background the entry is mapped to a SNOMED CT term. See


[http://help.visionhealth.co.uk/DLM700/Consultation\\_Manager/Content/ConMgr/General/Vision\\_Local\\_Codes.htm](http://help.visionhealth.co.uk/DLM700/Consultation_Manager/Content/ConMgr/General/Vision_Local_Codes.htm) for details.

### QOF




-  **Vision+ England - Vision+ Reports and Templates** are updated to **QOF Version 44** for practices in England.
-  **Vision+ Wales** – For practices in Wales using **QOF Version 37**, the reference date has changed from April 1<sup>st</sup> to October 1<sup>st</sup>.
-  **QOF Influenza Wales** – The QOF rules relating to Influenza are updated in this latest release to accommodate the new date ranges.
-  **CHD005 England (Patients taking aspirin, clopidogrel or anti-coagulant)** – Previously, when exception coding patients for all three drug areas, the alert prompt and template remained as if it was still outstanding (black text). This is resolved with QOF Version 44.
-  **Hypertension HYP003 and HYP004 England (Patients with a blood pressure 140/90 and 150/90 or less)** – The Blood pressure button was missing from the template. This is resolved with QOF Version 44.
-  **Vision+ Contract Manager England** – The Cervical screening indicators CS005 and CS006 within **Contract Manager**, only showed

the maximum threshold required to meet the target on the report. From SIS 10500, the report now displays both the minimum and maximum threshold.


### Vision+ Calculators

-  **eGFR** – The eGFR graph is no longer available because it was using the wrong weight for Cockcroft and Gault.


### Vision+ Templates

-  **History tab (additional information)** – The history tab within a template now displays additional information. See [Vision+ Templates](#).
-  **History tab (date)** – When viewing medication on the history tab within a template, the date the prescription was last issued displays. However, if the prescription had not been issued the date displayed as 30/12/1899. From SIS 10500 the date column is now blank.
-  **Smoking** – Smoking was missing from some Vision+ templates. This is corrected in SIS 10500.


### Vision+ Reports

-  **electronic Frailty Index (eFI) Calculator** – Previously, when running the eFI bulk calculator add option from **Practice Reports**, an error displayed. This is resolved in SIS 10500.


### Vision+ Template Designer

-  **Controls Toolbar** – The **Controls** toolbar within **Template Designer** includes some extra tools which include current repeats, acutes, repeat usage and problems. See [Vision+ Template Design](#).

### Anticoagulation Manager

-  **Warfarin Monitoring Reports** – The **Patient Scan** option within the Warfarin monitoring reports now runs a lot faster when identifying patients who need a diagnosis indication.

### Outcomes Manager


-  **Rule Designer** – Writing a pathway to search for eGFR>59 did not include eGFR>59 when searching in Vision 3. This is resolved in SIS 10500.

# Vision+ Templates

All **Vision+ templates** include a **history** tab. The history tab enables you to view a variety of clinical information for a patient, not just that relevant to the template. From any Vision+ template, select the history tab and see a categorised view of different types of information.

The categories now include additional columns. Below are some examples:

1. From **Consultation Manager**, select a **Vision+ Template** from the floating toolbar.
2. Select the **History** tab.

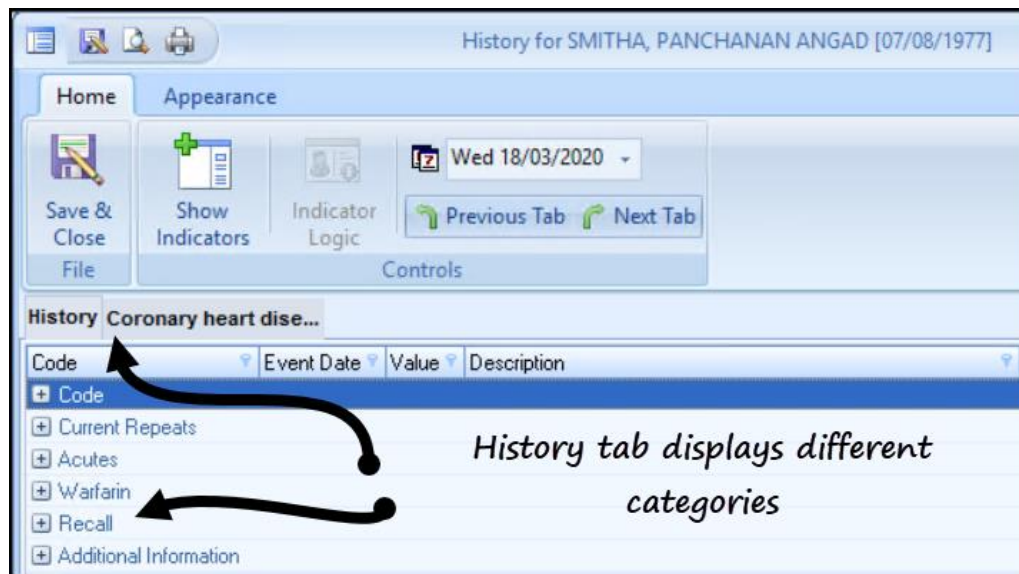

3. To view the data in the different categories select **expand** .


## Viewing data

Once a category is expanded scroll to view the following additional columns:

*For medication you can now see Dosage, Quantity and Pack*


Code	Event Date	Value	Description	Staff	Start	Comments	Priorit	Dosage	Quar	Pack	Ran
+ Code											
+ Current Repeats											
+ Acutes											
325974004	06/11/2019		Alendronic acid 10mg tablets AELL		06.06.02.00			1 TABLET ONCE A D 28		tablet	
29927500100002	28/12/2014		AGRIPPAL VACCINE prefiller AELL		00.00.00.00			20141228	0		
29927500100002	11/02/2014		BECLOMETASONE DIPROF AELL		00.00.00.00			20140211	0		
29927500100002	13/09/2013		PANDEMRIX INFLUENZA A AELL		00.00.00.00			20130913	0		
+ Warfarin											
+ Recall											
+ Additional Information											

Coded entries now display if they are recorded as a Problem (it also itemises if it is an active problem)

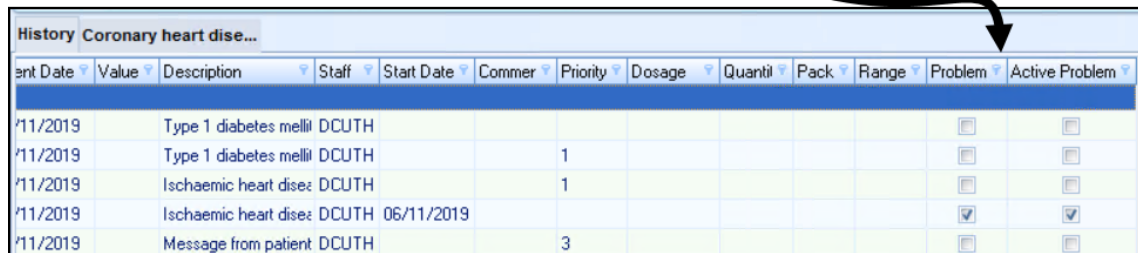



Event Date	Value	Description	Staff	Start Date	Commer	Priority	Dosage	Quantil	Pack	Range	Problem	Active Problem
11/2019		Type 1 diabetes mellit DCUTH									<input type="checkbox"/>	<input type="checkbox"/>
11/2019		Type 1 diabetes mellit DCUTH				1					<input type="checkbox"/>	<input type="checkbox"/>
11/2019		Ischaemic heart dise: DCUTH				1					<input type="checkbox"/>	<input type="checkbox"/>
11/2019		Ischaemic heart dise: DCUTH 06/11/2019									<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>
11/2019		Message from patient DCUTH				3					<input type="checkbox"/>	<input type="checkbox"/>

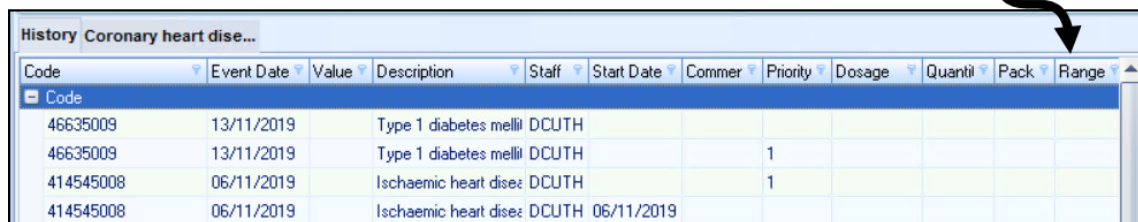Lab results where the > (greater than) option is used now display a range



Code	Event Date	Value	Description	Staff	Start Date	Commer	Priority	Dosage	Quantil	Pack	Range
46635009	13/11/2019		Type 1 diabetes mellit DCUTH								
46635009	13/11/2019		Type 1 diabetes mellit DCUTH				1				
414545008	06/11/2019		Ischaemic heart dise: DCUTH				1				
414545008	06/11/2019		Ischaemic heart dise: DCUTH 06/11/2019								


**Training Tip** – You can sort and filter patient information by

selecting the **filter**  option. To print what you are viewing select **print preview** or **print** .



# Vision+ Template Designer


Some additional controls have been added to **Template Designer**.

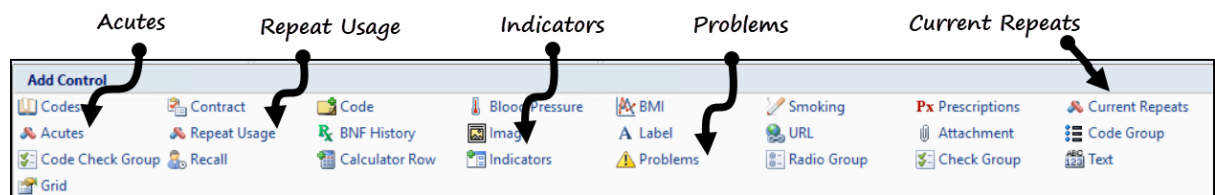


**Note** – Indicators are used in Rule Designer by the Pathway developers, and although visible in Template Designer, do not perform any function.


## Adding Active Repeats to a Template


You may want to view a patient's repeat medication within a template. This includes a history of the last ten repeat issues (however, this can be changed). Below is an example of how active repeats display in a Vision+ Template in Consultation Manager:

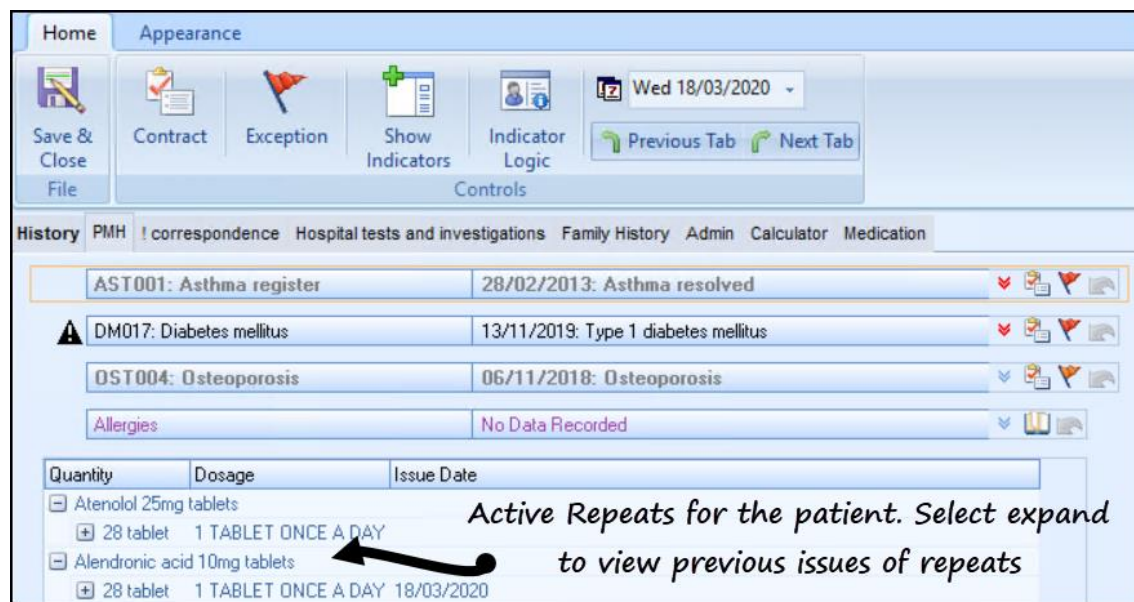



To add a **Current Repeats** row to a template:


1. From **Template Designer**, on the **Home** tab, select **Current Repeats** 
2. The **Current Repeats** row is added to the template:

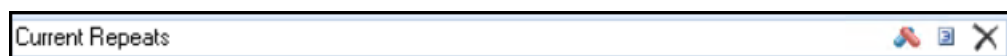


3. To limit the number of issues displayed in the template, select **Maximum number of rows** .
4. Select the number of rows to display and select **OK**.

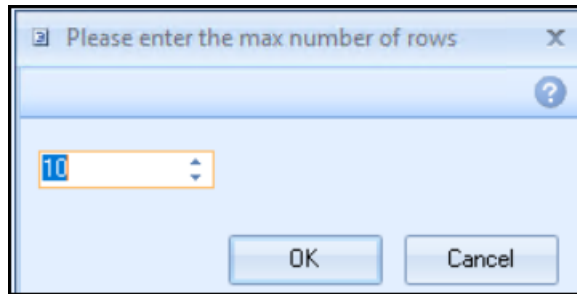

5. By default, only active repeats display in the template. To include inactive repeats select **Show inactive repeats**  at the end of the row.
6. **Show inactive repeats** displays. Select **Yes** and **OK** to view inactive repeats in the template.

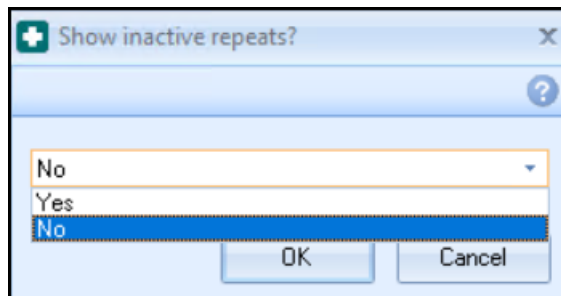

---

**Note** – The  button at the end of the row is used in Vision's Rule Designer module, and although visible does not perform any function in Template Designer.

---


## Adding an Acute Row

You can now view a patient's acute medication history within a template. By default, the last ten acute issues are listed, however, this can be changed. Below is an example of how acutes display in a Vision+ template:


The screenshot shows a patient's medical history template with tabs for History, PMH, correspondence, Hospital tests and investigations, Family History, Admin, Calculator, and Medication. The Medication tab is active, displaying a list of conditions: AST001: Asthma register (28/02/2013: Asthma resolved), DM017: Diabetes mellitus (13/11/2019: Type 1 diabetes mellitus), OST004: Osteoporosis (06/11/2018: Osteoporosis), and Allergies (No Data Recorded). Below this is a table of acute medications with columns for Quantity, Dosage, and Issue Date. The table lists Atenolol 25mg tablets (28 tablet, 1 TABLET ONCE A DAY) and Alendronic acid 10mg tablets (28 tablet, 1 TABLET ONCE A DAY, 18/03/2020). An arrow points to the 'Acute medication' section.

Quantity	Dosage	Issue Date
28 tablet	1 TABLET ONCE A DAY	
28 tablet	1 TABLET ONCE A DAY	18/03/2020


To add an **Acute** row to a template:

1. From **Template Designer**, on the **Home** tab, select **Acutes**  from the **Controls** area of the toolbar.
2. The **Acutes** row is added to the template:

The screenshot shows a single row in a template labeled 'Acutes' with a small icon at the end of the row.

3. To change the number of acutes displayed in the template, select **Maximum number of rows** .
4. Select the number of acute rows to display and select **OK**.

The screenshot shows a dialog box titled 'Please enter the max number of rows' with a spinner set to 10 and OK/Cancel buttons.


**Note** – The  button at the end of the row is used in Vision's Rule Designer module, and although visible does not perform any function in Template Designer.

## Adding a Repeat Usage Row


You can now view a patient's repeat usage in the last one year to assist monitoring compliance.

The following displays:


- For each repeat item, it looks to see what was prescribed in the past year.
- If the drug was prescribed earlier than one year ago, the bar is green. If it was prescribed sometime within the past year, it is a white bar up to the very first prescription that was issued, and then green from that point onwards.
- Each time a prescription is generated within that one year, it displays a black line.
- For a patient most of the bars will be green, with an equal spacing of black lines along the bar. This means the patient has been regularly collecting their prescription monthly, quarterly, etc. If the patient is collecting the prescription too frequently the black bars will be close together.

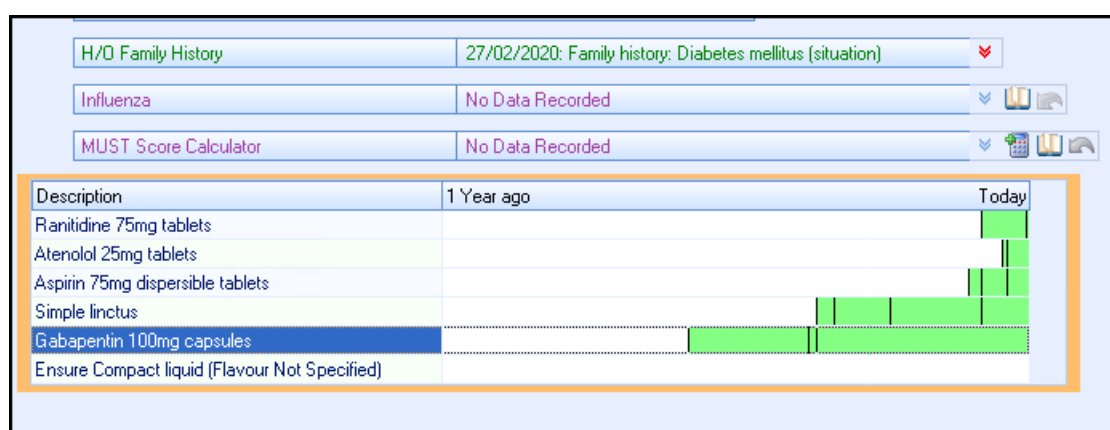

To add a **Repeat Usage** row to a template:


1. From **Template Designer**, on the **Home** tab, select **Repeat Usage** 
2. The **Repeat usage** row is added to the template:



3. To change the number of rows in the template, select **Maximum number of rows** .
4. Select the number of rows to display and select **OK**.

The screenshot shows a dialog box titled 'Please enter the max number of rows'. It has a text input field with the number '10' entered. There are 'OK' and 'Cancel' buttons at the bottom.

---

**Note** – The  button at the end of the row is used in Vision's Rule Designer module, and although visible does not perform any function in Template Designer.

---


## Adding a Problem Row

You can now view a patient's problems within a template. By default, a patient's active problems display but can be changed to include inactive problems.


Below is an example of how problems display in the template:


History	PMH	! correspondence	Hospital tests and investigations	Family History	Admin	Calculator	Medication
AST001: Asthma register		28/02/2013: Asthma resolved					
DM017: Diabetes mellitus		13/11/2019: Type 1 diabetes mellitus					
OST004: Osteoporosis		06/11/2018: Osteoporosis					
Description		Start Date					
Ischaemic heart disease		06/11/2019					


*Problems display with the start date. In Template design if inactive problems are selected, the end date also displays*


1. To add a **Problem** row to a template:
2. From **Template Designer**, on the **Home** tab, select **Problems**  from the **Controls** area of the toolbar.
3. The **Problem** row is added to the template:



3. To change the number of problem rows displayed in the template, select **Maximum number of rows** .
4. Select the number of rows to display and select **OK**.

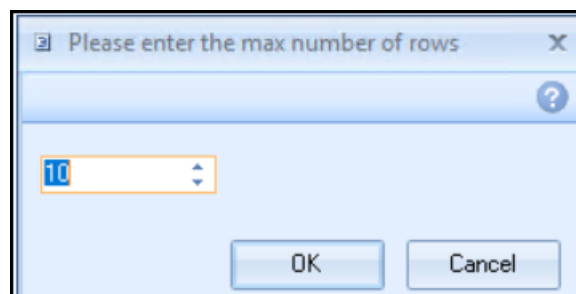

5. To view inactive problems within the template select **Show inactive problems** .
6. Select **Yes** and OK.

**Note** – The  button at the end of the row is used in Vision's Rule Designer module, and although visible does not perform any function in Template Designer.